




Medscope
Biotech

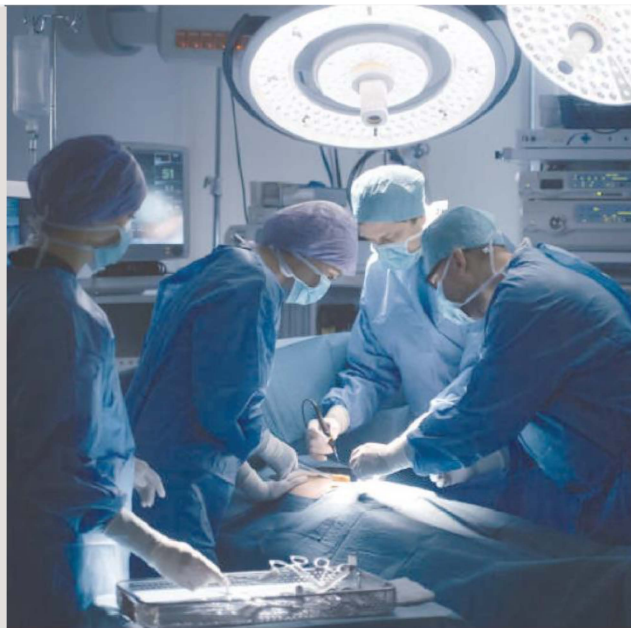

The Best Choice

Medscope Minimally Invasive Surgery Medical Device




SOLID CLIP™

CA16 Series 10mm

Single Use Clip Applier



SOLID CLIP™ applier is a disposable single-use sterile surgical instrument which provides the modern minimally invasive surgery a best occlusion solution

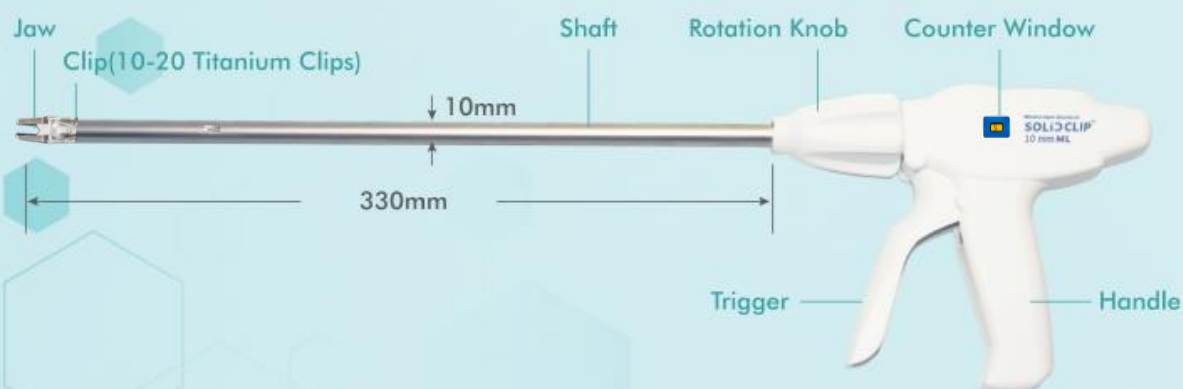

Specification

Model No.	Shaft Diameter	Clip Size	Number of Clips	Packaging
CA165	10mm	Large	10	5pcs/box, 4 boxes/carton
CA166	10mm	Large	15	5pcs/box, 4 boxes/carton
CA167	10mm	Medium/Large	15	5pcs/box, 4 boxes/carton
CA168	10mm	Medium/Large	20	5pcs/box, 4 boxes/carton
CA169	10mm	Medium/Large	10	5pcs/box, 4 boxes/carton

Clip Dimensions

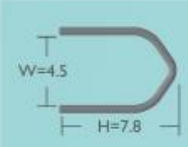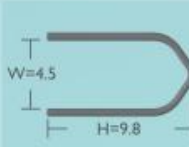Dimension Unit : mm	Clip Size		Medium / Large Size (ML)		Large Size (L)	
	(ML)	(L)	Unit : mm		Unit : mm	
W	4.5	4.5				
H	7.8	9.8				
L	9	11				

Features

Jaw Design



- 15° slop design for better visibility

Jaw Design



- Anti-drop design of the Jaw

Clip Design



- Inner-crossed pattern of the clip offers a firm hold on the vessel

Clip Design



- The closing feature of the clip can provide a proper occlusion of vessel

Clip Material



- Implantable CP Gr.1 titanium clips.
The compositions is ASTM F67 Standard Compliant

Clip Counter Window



- Preload 10~20 titanium clips
- The yellow block shown means that there are 5 clips left in this product. "0" means that there are ZERO clips left and the trigger will stop functioning.

Rotation Knob



The rotation knob can make the shaft rotating 360° in either direction

Handle



Ergonomic handle design brings excellent control from precise clip loading to closure.


Trigger



- Anti-slip design of the trigger
- Anti-fire design of the trigger when running out of the clips

SOLiDCLIP™ 5 mm

Single-use Clip Applier

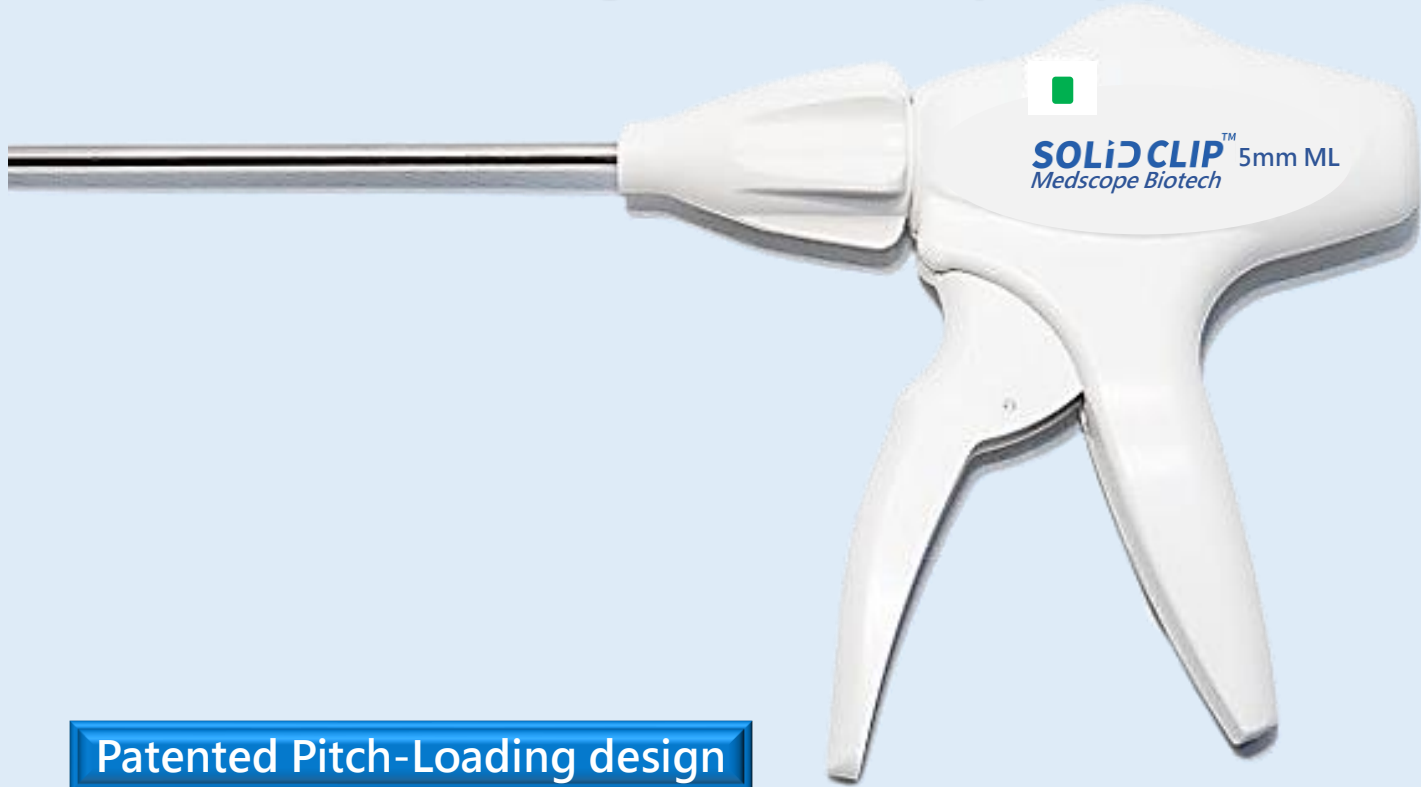

Patented Pitch-Loading design



Anti-drop Design

of the
Jaw Structure.




Pitch-Loading
Clip feeding Design

Lighter & Stable
Ligating Force.



Interlock Design

when running out
of the clips.

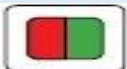Model					Preloaded Clip Number	Color Indicator			Dimensions
	Clip Size	W	mm H	L		Clip Number			
									
CA-568	M/L	3.5	8.2	9.0	15 Clips	1st ~ 12th	13th ~ 15th	No Clips	5 x 330 mm

Medscope Biotech Co., Ltd.

2F, No. 8, Keyi Street, Miaoli County, 350 Taiwan, R.O.C.

Tel: +886-37-580055 www.medscope.com.tw, Email: globalsales@medscope.com.tw

SOLiD CLIP™



Medscope Laparoscopy



Medscope Biotech Co., Ltd.

2F, No. 8, Keyi Street, GuangYuan Technology Park, Zhunan Town, 35059 Miaoli County, Taiwan

Tel: +886 (037) 580055 Fax: +886 (037) 585050

E-mail: globalsales@medscope.com.tw

www.medscope.com.tw